


RE-TUBING YOUR EARMOULD

1. Detach the hearing aid from earmould by holding the elbow of the hearing aid and pulling gently on the earmould. Remove the old tubing from the earmould (put to one side as you will need it later).



2. Take a piece of tubing and cut the end into a point so you can thread this through the hole in the earmould.
3. Thread the pointed end of the tubing through from the outside of the earmould. Pull until the bend in the tube meets the earmould and the tubing "stands up".

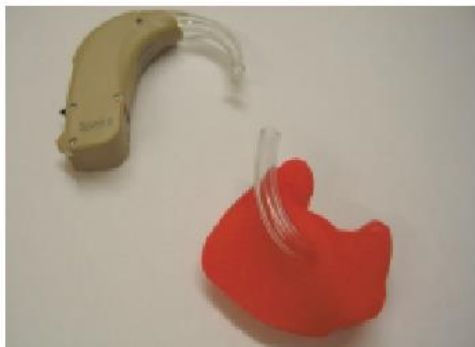

4. Cut the excess tubing off at the tip of the earmould. Line up the old piece of tubing against the earmould and use this as a guide to cut the other end of the tubing to the correct length.



5. Now your new tube should look like this. Reattach the earmould to the hearing aid.